

◆ **EC 300**
ELECTROSURGERY





Mains Voltage: 220V AC \pm 10%, 50 Hz (110 V AC Optional)

Maximum Power Consumption: 800 V.A

Protection Class: Class I

Type of output: CF

Weight: 8.8Kg

Dimensions (W.H.D): 37 x 16 x 46 cm


EC 300 offers various operating modes for open surgery, as well as laparoscopy. High power for cutting initiation in TUR achieves excellent results in fluid medium in urology or arthroscopy.

EC 300 is an electrosurgical unit for general surgery, Bipolar coagulation and Bipolar cutting in saline solution.

It also offers standard monopolar cutting and coagulation modes which are necessary in conventional procedures.

FEATURES:

- Bipolar Cut
- Bipolar Coagulation (Auto-start Bipolar & Manual Bipolar)
- Monopolar Cut (Pure)
- Monopolar Cut (Blend1)
- Monopolar Cut (Blend2)
- Monopolar Cut (Blend3)
- Monopolar TUR Function
- Monopolar Coagulation (Swift)
- Monopolar Coagulation (Forced)
- Monopolar Coagulation (Spray)
- Monopolar Coagulation (Soft)
- Bipolar TUR Mode (Optional)

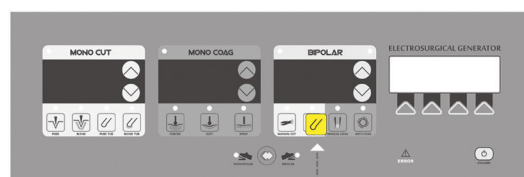

At Bipolar TUR, The Bipolar output power can be increased up to 250 Watts. So Cutting in Normal Saline is available



EC300
ELECTROSURGERY

